

PVAD Pulsatile VAD Tester

PVAD™ Pulsatile VAD Tester

The patented¹ PVAD™, a pulsatile ventricular assist device tester, is a real time tester that closely simulates all the physical conditions necessary to fully test artificial hearts and assist devices.

Overview


With the use of compliant aortic sections, the fluid dynamics of this tester can emulate *in vivo* conditions for flow rates and pressures. The test fixture provides mechanical equivalents for the left ventricle, aortic valve, compliant aorta, atrium, and mitral valve, completing a closed cardiovascular loop. A piston attached to a linear actuator with a servomotor provides the pumping function. This tester can be adjusted to test ventricular

assist devices or total artificial hearts by adjusting the volume per pulse and overall pressures and flows of the system. The sample chamber is temperature-controlled and allows the devices to be exposed to saline at body temperature for the duration of the test. The PVAD™ controlling software is written in the LabVIEW® graphical language which allows for simultaneously monitoring several conditions, managing operating parameters at programmable intervals, and collecting data from the instrumentation. Multiple testers can be linked together and operated from a central source.

1 US Patent # 7,621,192 B2

PVAD Tester Specifications

- Overall Dimensions: 18"L x 18"W x 71"H (minus control system)
- Power Requirements: 115V AC
- Maximum Fluid Displacement: 300mL
- Number of Samples: up to 3
- Volume Capacity: 28.4 liters
- Temperature Range: ambient to 40°C
- Accuracy of Temperature: -1°C +0.3°C
- Speed: Computer controlled stroke from 1 to 120 bpm
- Dry Weight: 203lbs

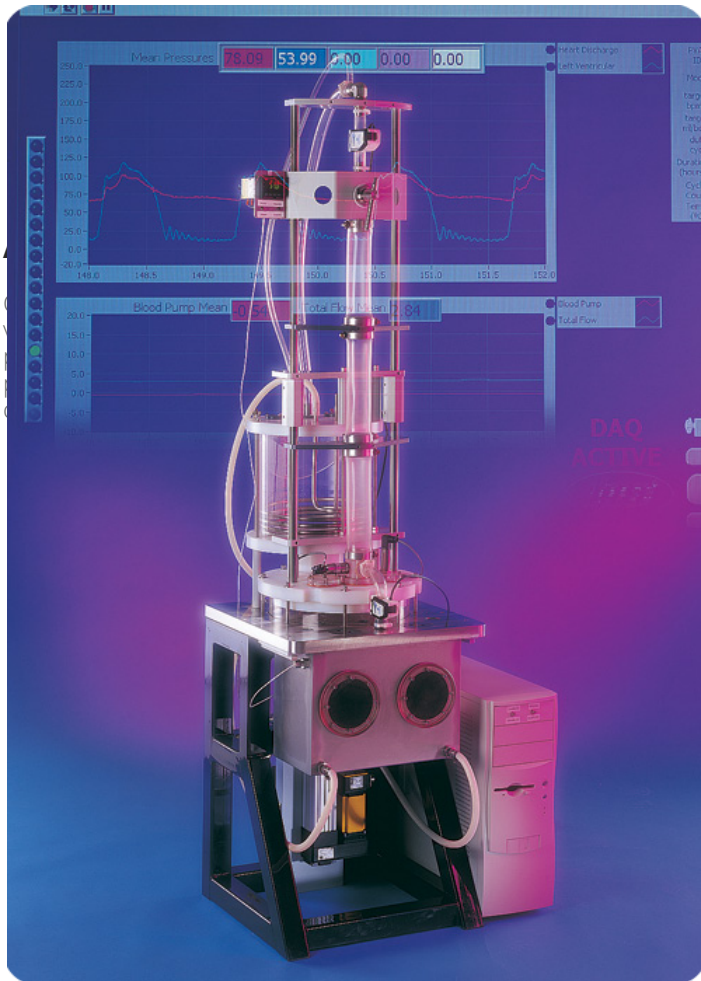

Advantages of the PVAD Tester

- Programmable for several different physiological states.
- Temperature controlled saline bath for samples.
- Can test up to three samples.
- Safety switches for temperature and pressure.
- Compliant aorta to allow for setting physiological flows.

To receive a customized proposal, contact us today at:

SD 534
111516

Dynatek Labs, Inc.

105 East 4th Street
Galena, MO 65656
800.325.8252
1.417.357.6155
www.dynateklabs.com
salesdd@dynateklabs.com

